


Installation Instructions

for the LRLIGHT POST

PROCEDURE:

1. Remove the mounting bracket cover by loosening the four cover screws (Refer to Illustration).
2. Position the light post assembly on the mounting bracket and screw the 3/4" bolt and washer into the tapped hole. Tighten finger tight leaving the light post assembly loose enough to allow the leveling screws to be adjusted.
3. Using a spirit level against the top portion of the light post, loosen or tighten the leveling screws to level the post. Be sure to check the level in all directions by moving the spirit level in all directions around the post.
4. Tighten the 3/4" bolt and check the level again. Repeat the leveling procedure if necessary.

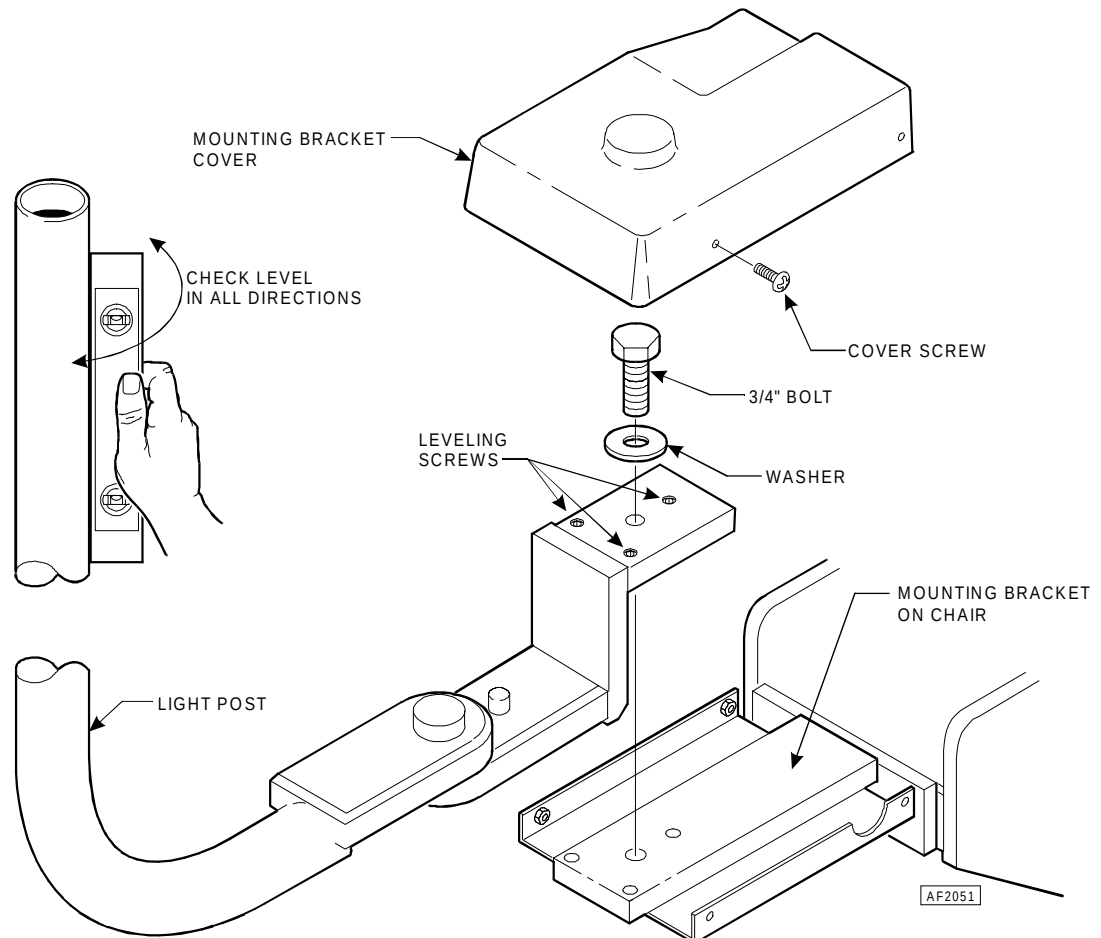

Knight®
by Midmark

Midmark Corporation • P O Box 286 • 60 Vista Drive
Versailles • Ohio • 45380-0286 • www.midmark.com
Phone: 937-526-3662 Fax: 937-526-5542

TPW1783 Rev. B